# NADAL<sup>®</sup> COVID-19 Ag Test (test cassette)

REF 243103N-25PA

<b>de</b>	Gebrauchsanweisung	2	<b>fi</b>	Käyttöohje	42
<b>en</b>	Instructions for use	7	<b>sv</b>	Användarinstruktioner	47
<b>fr</b>	Instructions d'utilisation	12	<b>da</b>	Brugervejledning	52
<b>es</b>	Instrucciones de uso	17	<b>no</b>	Bruksanvisning	57
<b>it</b>	Istruzioni per l'uso	22	<b>nl</b>	Gebruiksaanwijzing	62
<b>pl</b>	Sposób użycia	27	<b>ro</b>	Instrucțiuni de utilizare	67
<b>pt</b>	Instruções de Utilização	32	<b>bg</b>	Инструкции за употреба	72
<b>cs</b>	Návod k použití	27		Symbols	78
				Our Teams	80

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#### 15. Referenzen

- Cui J, U F, Shi ZL, Origin and evolution of pathogenic coronaviruses, Nat Rev Microbiol 2019; 17:181-192.

### 1. Intended Use

The NADAL® COVID-19 Ag Test is a lateral flow chromatographic immunoassay for the qualitative detection of SARS-CoV-2 viral nucleoprotein antigens in human nasal, nasopharyngeal or oropharyngeal specimens (see section 12 'Limitations'). This test is intended for use as an aid in the diagnosis of infections with SARS-CoV-2. Note that the concentration of viral nucleoprotein antigens may vary in the course of the disease and might fall below the detection limit of the test. Possible infectiousness of test subjects cannot be ruled out based on negative test results. The test procedure is not automated and requires no special training or qualification. The NADAL® COVID-19 Ag Test is designed for professional use only.

### 2. Introduction and Clinical Significance

COVID-19 (Coronavirus disease) is an infectious disease caused by the recently discovered coronavirus SARS-CoV-2. The most common symptoms of COVID-19 are a fever, dry cough, fatigue, sputum production, shortness of breath, a sore throat and headaches. Some patients may experience myalgia, chills, nausea, nasal congestion and diarrhoea. These symptoms begin gradually and are mild in most cases. Some individuals become infected but do not develop any symptoms and do not feel unwell. Most people (about 80%) recover from the disease without special treatment. Approximately one in six people who get infected with COVID-19 become seriously ill and develop difficulty breathing. Elderly people, and those with pre-existing conditions, such as high blood pressure, heart problems or diabetes, are more likely to develop serious illness.

COVID-19 is transmitted via respiratory droplets that are exhaled by infected people via coughing, sneezing or talking. These droplets can be inhaled or ingested directly by other people or can contaminate surfaces, which can then be infectious for several days. Most estimates of the incubation period for COVID-19 range from 1 to 14 days, during which people might already be infectious without showing disease symptoms.

### 3. Test Principle

The NADAL® COVID-19 Ag Test is a lateral flow chromatographic immunoassay for the qualitative detection of SARS-CoV-2 viral nucleoprotein antigens in human nasal, nasopharyngeal or oropharyngeal specimens.

Anti-SARS-CoV-2 antibodies are immobilised in the test line region (T) of the membrane. A specimen is added to an extraction tube containing buffer in order to release SARS-CoV-2 antigens. During the test, extracted antigens bind to anti-SARS-CoV-2 antibodies conjugated to coloured particles and pre-coated onto the sample pad of the test cassette. The mixture then migrates along the membrane chromatographically by capillary action and interacts with the reagents on the membrane. The complexes are then captured by anti-SARS-CoV-2 antibodies in the test line region (T). Excess coloured particles are captured in the control line region (C). The presence of a coloured line in the test line region (T) indicates a positive result. The absence of a coloured line in the test line region (T) indicates a negative result.

The formation of a coloured line in the control line region (C) serves as a procedural control, indicating that the proper volume of specimen has been added and membrane wicking has occurred.

### 4. Reagents and Materials Supplied

#### • 25 NADAL® COVID-19 Ag test cassettes\*

- Additional material provided according to 93/42/EEC: Due to possible supply shortages of COVID-19 related accessory medical products, the swab manufacturer might change. Therefore, the supplied swabs are from one of the manufacturers listed below.

#### a) 25 sterile swabs, CE 2797

Puritan Medical Products Company LLC  
31 School Street

Guilford, Maine 04443-0149 USA (authorised EU representative EMERGO EUROPE, The Hague, The Netherlands)

#### b) 25 sterile swabs, CE 0197

Jiangsu Changfeng Medical Industry Co., Ltd  
Touqiao Town, Guangling District, Yangzhou,  
Jiangsu 225109 China (authorised EU representative Lins Service & Consulting GmbH, Obere Seegasse 34/2, 69124 Heidelberg, Germany)

#### c) 25 sterile swabs, CE 0197

CITOTEST LABWARE MANUFACTURING CO., LTD  
No.48, Xinxu Road, Haimen, Jiangsu province  
(authorised EU representative WellKang Ltd, Enterprise Hub, NW Business Complex, 1 Beraghmore Rd., Derry, BT48 8SE, Northern Ireland)

#### d) 25 sterile swabs, Copan Floqswabs; CE 0123

Copan Italia S.p.A.,  
Via Perotti 10, 25125 Brescia, Italy

- 25 extraction tubes containing buffer (400 µL each)\*\*, incl. dropper caps
- 1 reagent holder
- 1 package insert

\*containing the preservative sodium azide: <0.1%

\*\*Buffer containing the following preservative: sodium azide: <0.1 mg/mL

The detergents contained in the buffer lyse and neutralise the virus.

No hazard labelling is required according to Regulation (EC) N° 1272/2008 CLP. Concentrations are below exemption threshold.

### 5. Additional Materials Required

- Timer

### 6. Storage & Stability

Test kits should be stored at 2-30°C until the indicated expiry date. Test cassettes are stable until the expiry date printed on the foil pouches. Test cassettes must remain in the sealed foil pouches until use. Do not freeze the test kit. Do not use tests beyond the expiry date indicated on the packaging. Care should be taken to protect test kit components from

contamination. Do not use test kit components if there is evidence of microbial contamination or precipitation. Biological contamination of dispensing equipment, containers or reagents can lead to inaccurate results.

## 7. Warnings and Precautions

- For professional *in-vitro* diagnostic use only.
- Carefully read through the test procedure prior to testing.
- Do not use the test beyond the expiration date indicated on the packaging.
- Do not use test kit components if the primary packaging is damaged.
- Tests are for single use only.
- Do not add specimens to the reaction area (result area).
- In order to avoid contamination, do not touch the reaction area (result area).
- Avoid cross-contamination of specimens by using a new extraction tube for each specimen obtained.
- Do not substitute or mix components from different test kits.
- Do not use the buffer if it is discoloured or turbid. Discolouration or turbidity may be a sign of microbial contamination.
- Do not eat, drink or smoke in the area where specimens and test kits are handled.
- Wear appropriate personal protective equipment, such as face mask, isolation gown, gloves and eye protection during specimen collection, preparation and test procedure.
- Handle all specimens as if they contain infectious agents. Observe established precautions for microbiological risks throughout all procedures and standard guidelines for the appropriate disposal of specimens.
- Further specimen processing and patient management should follow local COVID-19 guidelines and regulations.
- The test kit contains products of animal origin. Certified knowledge of the origin and/or sanitary state of the animals does not completely guarantee the absence of transmissible pathogenic agents. It is therefore recommended that these products be treated as potentially infectious and handled in accordance with usual safety precautions (e.g., do not ingest or inhale).
- Temperature can adversely affect test results.
- Used testing materials should be disposed of according to local regulations.

## 8. Specimen Collection and Preparation

### Nasal specimen:

- It is important to obtain as much secretion as possible. Insert the swab into the nostril.
- Gently push the swab until resistance is met at the level of the turbinate (up to 2.5 cm into the nostril).
- Rotate the swab 5 times against nasal wall to ensure that both mucus and cells are collected.
- Slowly withdraw the swab while continuing to rotate it.
- Repeat this process for the other nostril using the same swab in order to ensure that sufficient specimen volume is collected from both nasal cavities.

### Oropharyngeal specimen:

- Gently insert a sterile swab into the pharynx and collect secretions by brushing the swab several times against the reddened posterior pharyngeal wall and both tonsillar pillars. Avoid touching the tongue, teeth and gums.

### Nasopharyngeal specimen:

- Insert the swab into the nostril, parallel to the palate (not upwards) until resistance is encountered or the distance is equivalent to that from the ear to the nostril of the patient, indicating contact with the nasopharynx.
- Gently rub and roll the swab. Leave swab in place for several seconds to absorb secretions.
- Slowly remove the swab while rotating it. Specimens can be collected from both nostrils using the same swab, but it is not necessary to collect specimens from both sides if the tip is saturated with fluid from the first collection.


### Note:

- Use only synthetic fibre swabs with plastic shafts. Do not use calcium alginate swabs or swabs with wooden shafts, as they may contain substances that inactivate some viruses and inhibit further testing.
- Swab specimens should be tested immediately after collection. Use freshly collected specimens for best test performance.
- If not tested immediately, swab specimens can be stored in viral transport media (VTM) without denaturing agents at 2-8°C for 24 hours after collection.
- Viral transport media (VTM) without denaturing agents as well as VTM containing virus-inactivating agents and evaluated with the test can be used during specimen storage for the subsequent antigen detection using the NADAL® COVID-19 Ag Test. In order to influence the sensitivity as little as possible, a low volume of VTM (max. 1 mL) is recommended. The test procedure when using viral transport media (VTM) is available separately:
  - under the following link:  
<https://www.nal-vonminden.com/en/vtm-ifu>
  - by request at +49 941 29010-0 or  
[info@nal-vonminden.com](mailto:info@nal-vonminden.com)
- Do not use specimens that are obviously contaminated with blood, as it may interfere with the flow of specimens and lead to inaccurate test results.

## 9. Test Procedure

**Bring tests, specimens, extraction tubes containing buffer and/or controls to room temperature (15-30°C) prior to testing.**


1. Hold an extraction tube containing buffer and labelled with the patient or control identification vertically. Peel off the aluminum foil from the tube.
2. Insert the swab with the collected specimen into the tube. Swirl the swab and squeeze it at least 10-15 times by compressing the wall of the extraction tube against the swab to extract the antigens contained in the swab.
3. Remove the swab, pressing it firmly against the wall of the tube to release

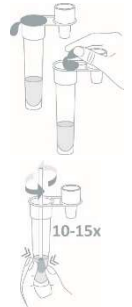as much liquid as possible. Dispose of the swab in accordance with guidelines for the handling of infectious agents.

- Attach a dropper cap to the tube.
- Place the tube into the designated area of the reagent holder.
- Remove the test cassette from the foil pouch and use it as soon as possible. The best results will be obtained if the test is performed immediately after opening the foil pouch. Label the test cassette with the patient or control identification.
- Place the test cassette on a clean and level surface.
- Invert the tube and transfer 3 drops of the extracted solution to the specimen well (S) of the test cassette.
- Start the timer.
- Wait for the coloured line(s) to appear. Read the test result after 15 minutes. Do not interpret the result after more than 20 minutes.

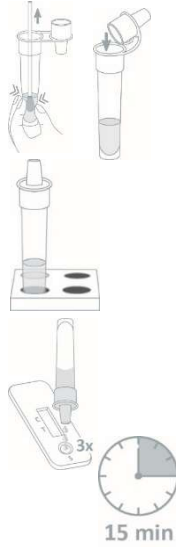A coloured line appearing in the control line region (C) is considered an internal procedural control. It confirms sufficient specimen volume, adequate membrane wicking and correct procedural technique.

*Good laboratory practice (GLP)* recommends the use of external control materials to ensure proper test kit performance.

## 12. Limitations

- The NADAL® COVID-19 Ag Test is for professional *in-vitro* diagnostic use only. It should be used for the qualitative detection of SARS-CoV-2 viral nucleoprotein antigens in human nasal, nasopharyngeal or oropharyngeal specimens only. Neither the quantitative value nor the rate of increase/decrease in the concentration of SARS-CoV-2 viral nucleoprotein antigens can be determined with this qualitative test.
- The NADAL® COVID-19 Ag Test only detects the presence of SARS-CoV-2 viral nucleoprotein antigens in specimens and should not be used as the sole criterion for a diagnosis of COVID-19.
- Both viable and non-viable SARS-CoV-2 viruses can be detected using the NADAL® COVID-19 Ag Test.
- The sections 'Specimen Collection and Preparation' as well as 'Test Procedure' must be followed closely while testing. Failure to follow them may lead to inaccurate test results because the antigen concentration in the swab is highly dependent on the correct procedure.
- As with all diagnostic tests, all results should be interpreted in conjunction with other clinical information available to the physician.
- In the course of SARS-CoV-2 infection, the concentration of viral nucleoprotein antigens may fall below the detection limit of the test.
- If the test result is negative and clinical symptoms persist, additional testing using other clinical methods is recommended. A negative result does not at any time preclude the possibility of a SARS-CoV-2 infection and should be confirmed via molecular assay.
- Positive and negative predictive values are highly dependent on prevalence. The local prevalence should be taken into consideration when interpreting diagnostic test results.
- Positive results do not preclude co-infections with other pathogens (e.g. influenza virus A/B).

## 13. Expected Values

SARS-CoV-2 viral particles are normally present in the respiratory tracts of COVID-19 patients. A positive test result can indicate an acute infection. Virus concentrations in nasal, nasopharyngeal or oropharyngeal swab specimens may vary in the course of the disease and might fall below the detection limit of rapid tests, even though patients are still showing symptoms. Conversely, the virus might continue to be detectable over long periods of time even in convalescent patients. Possible infectiousness of test subjects cannot be ruled out based on negative test results.


## 14. Performance Characteristics

**Clinical performance (specimens extracted in buffer without VTM)**

## 10. Result Interpretation

### Positive:

Two coloured lines appear in the result area. One line appears in the control line region (C) and the other line appears in the test line region (T).

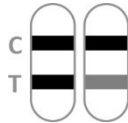**Note:** The colour intensity in the test line region (T) may vary depending on the concentration of SARS-CoV-2 viral nucleoprotein antigens in the specimen. Any shade of colour in the test line region (T) should be considered a positive result. Note that this is a qualitative test only and it cannot determine the analyte concentration in the specimen.

### Negative:

Only one coloured line appears in the control line region (C). No coloured line appears in the test line region (T).

### Invalid

The control line (C) fails to appear. Results from any test which has not produced a control line at the specified reading time must be discarded. Please review the procedure and repeat the test with a new test cassette. If the problem persists, discontinue using the test kit immediately and contact your distributor.

Insufficient specimen volume, incorrect operating procedure or expired tests are the most likely reasons for control line failure.

## 11. Quality Control

An internal procedural control is included in the test cassette:

### Diagnostic sensitivity and specificity

The NADAL® COVID-19 Ag Test was evaluated with clinical nasopharyngeal and oropharyngeal swab specimens whose status was confirmed using RT-PCR ( $C_t$  range positive: 20-37). The sensitivity was calculated for the range from high to medium viral load ( $C_t$  20-30) and from high to very low viral load ( $C_t$  20-37). The results are presented in the following tables.

		RT-PCR, $C_t$ 20-30		
		Positive	Negative	Total
NADAL® COVID-19 Ag Test	Positive	120	0	120
	Negative	3	726	729
	Total	123	726	849

Diagnostic sensitivity ( $C_t$  20-30): 97.6% (93.1% - 99.2%)\*

Overall agreement ( $C_t$  20-30): 99.7% (99.0% - 99.9%)\*

Diagnostic specificity: >99.9% (99.5% - 100%)\*

\*95% confidence interval

		RT-PCR, $C_t$ 20-37		
		Positive	Negative	Total
NADAL® COVID-19 Ag Test	Positive	150	0	150
	Negative	37	726	763
	Total	187	726	913

Diagnostic sensitivity ( $C_t$  20-37): 80.2% (73.9% - 85.3%)\*

Overall agreement ( $C_t$  20-37): 96.0% (94.5% - 97.1%)\*

Diagnostic specificity: >99.9% (99.5% - 100%)\*

\*95% confidence interval

In order to show the strong dependence of the diagnostic sensitivity on the viral load, the following table demonstrates the sensitivity for different  $C_t$  value ranges of the reference PCR:

$C_t$ range	Sensitivity
20 – 25	97.12%
20 – 30	97.56%
20 – 32	96.21%
20 – 35	85.71%
20 – 37	80.21%

Please note that  $C_t$  values may vary between different PCR systems at the same virus concentration.

In an additional study, the NADAL® COVID-19 Ag Test was evaluated with clinical nasal swab specimens whose status was confirmed using RT-PCR. The results are presented in the following table.

		RT-PCR, $C_t$ <30		
		Positive	Negative	Total
NADAL® COVID-19 Ag Test	Positive	96	0	96
	Negative	6	138	144
	Total	102	138	240

Diagnostic sensitivity ( $C_t$  <30): 94.12% (87.64% - 97.81%)\*

Overall agreement ( $C_t$  <30): 97.50% (94.64% - 99.08%)\*

Diagnostic specificity: >99.9% (97.36% - 100%)\*

\*95% confidence interval

The specificity of the NADAL® COVID-19 Ag Test was also evaluated with 105 nasal swab specimens from asymptomatic

test subjects. The NADAL® COVID-19 Ag Test identified 105 out of 105 specimens as true negative (specificity >99.9%).

### Detection limit

The detection limit of the NADAL® COVID-19 Ag Test is  $2 \times 10^{2.4}$  TCID<sub>50</sub>/mL and was determined with a SARS-CoV-2 control with a known virus titre.

The detection limit of the NADAL® COVID-19 Ag Test is 0.4 ng/mL for recombinant SARS-CoV-2 nucleoprotein.

### Interfering substances

The following substances, normally present in respiratory specimens or artificially introduced into the respiratory tract, were evaluated at the concentrations listed below and showed no interference with the NADAL® COVID-19 Ag Test.

Substance	Concentration	Substance	Concentration
3 OTC nasal sprays	10%	Guaiaicol glyceryl ether	20 mg/mL
3 OTC mouthwashes	10%	Mucin	1%
3 OTC sore throat liquids	10%	Mupirocin	250 µg/mL
4-acetamidophenol	10 mg/mL	Oxymetazoline	10 mg/mL
Acetylsalicylic acid	20 mg/mL	Phenylephrine	10 mg/mL
Albuterol	20 mg/mL	Phenylpropanolamine	20 mg/mL
Chlorpheniramine	5 mg/mL	Relenza® (zanamivir)	20 mg/mL
Dexamethasone	5 mg/mL	Rimantadine	500 ng/mL
Dextromethorphan	10 mg/mL	Tamiflu® (oseltamivir)	100 mg/mL
Diphenhydramine	5 mg/mL	Tobramycin	40 mg/mL
Doxylamine succinate	1 mg/mL	Triamcinolone	14 mg/mL
Flunisolide	3 mg/mL		

### Cross-reactivity and microbial/viral interference

SARS-CoV-2 negative and positive specimens spiked with the following pathogens were tested using the NADAL® COVID-19 Ag Test:

HCoV-HKU1, HCoV-OC43, HCoV-NL63, HCoV-229E, measles virus, *Streptococcus pneumoniae*, Epstein-Barr virus, *Bordetella parapertussis*, *Bordetella pertussis*, influenza A (H1N1) pdm09, influenza A (H3N2), influenza A (H5N1), influenza A (H7N9), influenza A (H7N7), influenza B Victoria lineage, *Haemophilus influenzae*, *Candida albicans*, *Mycobacterium tuberculosis*, respiratory syncytial virus, adenovirus, parainfluenza virus type 1, 2, 3, 4, human metapneumovirus, rhinovirus, coxsackievirus type A16, norovirus, mumps virus, *Legionella pneumophila*, *Mycoplasma pneumoniae*, *Chlamydia pneumoniae*, *Streptococcus pyogenes*, *Streptococcus agalactiae*, *Streptococcus* group C, *Staphylococcus aureus*.

No cross-reactivity or microbial/viral interference with the specimens was observed when tested using the NADAL® COVID-19 Ag Test.

**Precision****Repeatability and reproducibility**

Precision was established by testing 10 replicates of negative, low positive and high positive controls.

Reproducibility was established by testing triplicates of negative, low and high positive controls. Testing was performed by 3 operators using 3 independent NADAL® COVID-19 Ag test lots at 3 different sites on 5 separate days.

The NADAL® COVID-19 Ag Test demonstrated acceptable repeatability and reproducibility. The negative and positive values were correctly identified >99% of the time.

**15. References**

1. Cui J, Li F, Shi ZL, Origin and evolution of pathogenic coronaviruses, *Nat Rev Microbiol* 2019; 17:181-192.
2. Su S, Wong G, Shi W, et al, Epidemiology, genetic recombination, and pathogenesis of coronaviruses, *Trends Microbiol* 2016;24:490-502.
3. Weiss SR, Leibowitz JL, Coronavirus pathogenesis, *Adv Virus Res* 2011; 81:85-164.

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- a) 25 tamponi sterili, CE 2797

Puritan Medical Products Company LLC  
31 School Street  
Guilford, Maine 04443-0149 USA (rappresentante autorizzato EU EMERGO EUROPE, The Hague, The Netherlands)

- b) 25 tamponi sterili, CE 0197

Jiangsu Changfeng Medical Industry Co., Ltd  
Touqiao Town, Guangling District, Yangzhou,  
Jiangsu 225109 China (rappresentante autorizzato EU Lins Service & Consulting GmbH, Obere Seegasse 34/2, 69124 Heidelberg, Germany)

- c) 25 tamponi sterili, CE 0197

CITOTEST LABWARE MANUFACTURING CO., LTD  
No.48, Xinxu Road, Haimen, Jiangsu province  
(rappresentante autorizzato EU WellKang Ltd, Enterprise Hub, NW Business Complex, 1 Beraghmore Rd., Derry, BT48 8SE, Northern Ireland)

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a) 25 sterile vattinnene, CE 2797

Puritan Medical Products Company LLC  
31 School Street  
Guilford, Maine 04443-0149 USA (autorisert EU-representant EMERGO EUROPA, The Hague, The Netherlands)

b) 25 sterile vattinnene, CE 0197

Jiangsu Changfeng Medical Industry Co., Ltd  
Touqiao Town, Guangming District, Yangzhou,  
Jiangsu 225109 China (autorisert EU-representant Lins Service & Consulting GmbH, Obere Seegasse 34/2, 69124 Heidelberg, Germany)

c) 25 sterile vattinnene, CE 0197

CITOTEST LABWARE MANUFACTURING CO., LTD  
No.48, Xinxu Road, Haimen, Jiangsu provins (autorisert EU-representant WellKang Ltd, Enterprise Hub, NW Business Complex, 1 Beraghmore Rd., Derry, BT48 8SE, Northern Ireland)

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nal von minden GmbH